

## Anti-Ub antibody (40-89 aa) (STJ96158)

STJ96158

### GENERAL INFORMATION


<b>Product Type</b>	Primary antibodies
<b>Short Description</b>	Rabbit polyclonal antibody anti-Polyubiquitin-B and Polyubiquitin-C and Ubiquitin-40S ribosomal protein S27a and Ubiquitin-60S ribosomal protein L40 (40-89 aa) is suitable for use in Western Blot, Immunoprecipitation, Immunohistochemistry, Immunofluo
<b>Applications</b>	WB/IP/IHC/IF/ELISA
<b>Host/Source</b>	Rabbit
<b>Reactivity</b>	Human/Mouse/Rat

### PRODUCT PROPERTIES


<b>Clonality</b>	Polyclonal
<b>Clone ID</b>	
<b>Concentration</b>	1 mg/mL
<b>Conjugation</b>	Unconjugated
<b>Purification</b>	The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen.
<b>Dilution</b>	WB 1:500-1:2000
<b>Range</b>	IHC 1:100-1:300 ELISA 1:20000 IF 1:50-200
<b>Formulation</b>	Liquid in PBS containing 50% Glycerol, 0.5% BSA and 0.02% Sodium Azide.
<b>Isotype</b>	IgG
<b>Storage</b>	Store at -20°C for up to 1 year from the date of receipt, and avoid repeat freeze-thaw cycles.
<b>Instruction</b>	

### TARGET INFORMATION

<b>Gene ID</b>	6233 7311 7314 RPS27A UBA52 RS27A_HUMAN RL40_HUMAN UBB_HUMAN
<b>Immunogen</b>	The antiserum was produced against synthesized peptide derived from the human Ubiquitin at the amino acid range 40-89
<b>Immunogen</b>	40-89 aa
<b>Region</b>	
<b>Specificity</b>	Ub Polyclonal Antibody detects endogenous levels of Ub protein.
<b>Immunogen Sequence</b>	

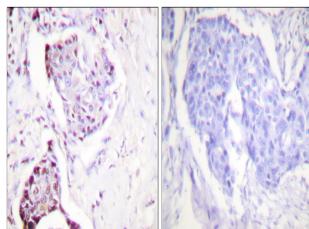

Immunohistochemistry analysis of paraffin-embedded human breast carcinoma tissue, using Ubiquitin Antibody. The picture on the right is blocked with the synthesized peptide.

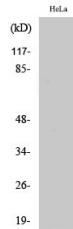

Western blot analysis of various cells using Ub Polyclonal Antibody. Secondary antibody was diluted at 1:20000

This product is suitable for in-vitro studies under the RESEARCH USE ONLY [RUO] licence. This product must not be used as for diagnostic or other medical purposes.  
St John's Laboratory Ltd, Knowledge Dock Business Centre, University Way, London, E16 2RD | Tel: 0208 223 3081