

Anti-Phospho-VIM-S39 antibody (STJ117904)

STJ117904

GENERAL INFORMATION


Product Type	Primary antibodies
Short Description	Rabbit polyclonal antibody anti-Phospho-VIM-S39 is suitable for use in Western Blot and Immunohistochemistry.
Applications	WB, IHC
Host/Source	Rabbit
Reactivity	Human, Mouse, Rat

PRODUCT PROPERTIES


Clonality	Polyclonal
Clone ID	
Concentration	
Conjugation	Unconjugated
Purification	Affinity purification
Dilution Range	WB 1:500-1:2000 IHC 1:50-1:100
Formulation	PBS containing 0.02% Sodium Azide, 50% Glycerol, pH7.3.
Isotype	IgG
Storage Instruction	Store in a freezer at -20°C and avoid freeze-thaw cycles.

TARGET INFORMATION


Gene ID	7431
Gene Symbol	VIM
Uniprot ID	VIME_HUMAN
Immunogen	A synthetic phosphorylated peptide around S39 of human VIM (NP_003371.2).
Immunogen Region	
Specificity	
Immunogen Sequence	

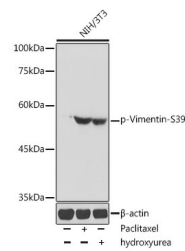

Western blot analysis of extracts of NIH/3T3 cells, using Phospho-Vimentin-S39 antibody (STJ117904) at 1:1000 dilution. NIH/3T3 cells were treated by Pacitaxel (100 nM/ml) at 37 °C for 20 hours. NIH/3T3 cells were treated by Hydroxyurea (4 mM) at 37 °C for 20 hours. Secondary antibody: HRP Goat Anti-rabbit IgG (H+L) at 1:10000 dilution. Lysates/proteins: 25ug per lane. Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Enhanced Kit. Exposure time: 1s.

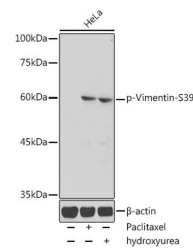

Western blot analysis of extracts of HeLa cells, using Phospho-Vimentin-S39 antibody (STJ117904) at 1:1000 dilution. HeLa cells were treated by Pacitaxel (100 nM/ml) at 37 °C for 20 hours. HeLa cells were treated by Hydroxyurea (4 mM) at 37 °C for 20 hours. Secondary antibody: HRP Goat Anti-rabbit IgG (H+L) at 1:10000 dilution. Lysates/proteins: 25ug per lane. Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Basic Kit. Exposure time: 180s.

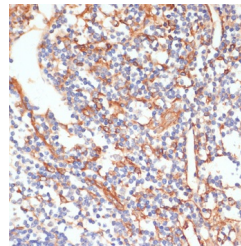

Immunohistochemistry of paraffin-embedded human appendix using Phospho-Vimentin-S39 antibody (STJ117904) at dilution of 1:100 (40x lens).

This product is suitable for in-vitro studies under the RESEARCH USE ONLY [RUO] licence. This product must not be used as for diagnostic or other medical purposes.
St John's Laboratory Ltd, Knowledge Dock Business Centre, University Way, London, E16 2RD | Tel: 0208 223 3081