

Anti-EFTUD2 antibody (1-205) (STJ29120)

STJ29120

GENERAL INFORMATION


Product Type	Primary antibodies
Short Description	Rabbit polyclonal antibody anti-EFTUD2 (1-205) is suitable for use in Western Blot, Immunohistochemistry, Immunofluorescence and Immunoprecipitation.
Applications	WB, IHC, IF, IP
Host/Source	Rabbit
Reactivity	Human, Mouse, Rat

PRODUCT PROPERTIES


Clonality	Polyclonal
Clone ID	
Concentration	
Conjugation	Unconjugated
Purification	Affinity purification
Dilution Range	WB 1:500-1:2000 IHC 1:50-1:200 IF 1:50-1:100 IP 1:20-1:50
Formulation	PBS containing 0.02% Sodium Azide, 50% Glycerol, pH7.3.
Isotype	IgG
Storage	Store in a freezer at -20°C and avoid freeze-thaw cycles.
Instruction	

TARGET INFORMATION


Gene ID	9343
Gene Symbol	EFTUD2
Uniprot ID	U5S1_HUMAN
Immunogen	Recombinant fusion protein containing a sequence corresponding to amino acids 1-205 of human EFTUD2 (NP_004238.3).
Immunogen Region	1-205
Specificity	
Immunogen Sequence	

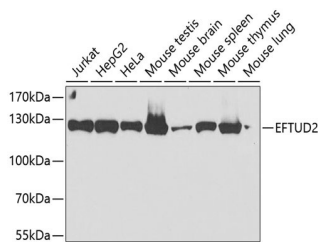


Western blot analysis of extracts of various cell lines, using EFTUD2 antibody (STJ29120) at 1:2000 dilution. Secondary antibody: HRP Goat Anti-rabbit IgG (H+L) at 1:10000 dilution. Lysates/proteins: 25ug per lane. Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Enhanced Kit. Exposure time: 20s.

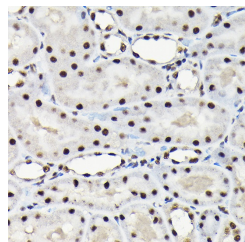

Immunohistochemistry of paraffin-embedded rat kidney using EFTUD2 rabbit polyclonal antibody (STJ29120) at dilution of 1:100 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6, 0 before commencing with immunohistochemistry staining protocol.

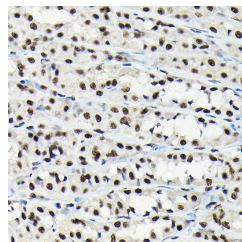

Immunohistochemistry of paraffin-embedded rat stomach using EFTUD2 rabbit polyclonal antibody (STJ29120) at dilution of 1:100 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6, 0 before commencing with immunohistochemistry staining protocol.

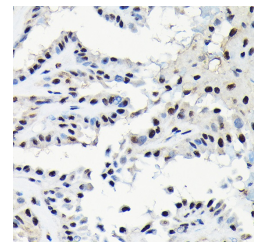

Immunohistochemistry of paraffin-embedded human thyroid cancer using EFTUD2 rabbit polyclonal antibody (STJ29120) at dilution of 1:100 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6, 0 before commencing with immunohistochemistry staining protocol.

This product is suitable for in-vitro studies under the RESEARCH USE ONLY [RUO] licence. This product must not be used as for diagnostic or other medical purposes.
St John's Laboratory Ltd, Knowledge Dock Business Centre, University Way, London, E16 2RD | Tel: 0208 223 3081