

Anti-Wiskott-Aldrich Syndrome/WASP antibody (Internal) (STJ70672)

GENERAL INFORMATION

Product Type	Primary antibodies
Applications	Pep-ELISA/IHC/IF/FC
Host / Source	Goat
Reactivity	Human/Mouse/Rat/Dog

PRODUCT PROPERTIES


Clonality	Polyclonal
Concentration	0.5 mg/mL
Conjugation	Unconjugated
Purification	Purified from goat serum by ammonium sulphate precipitation followed by antigen affinity chromatography using the immunizing peptide.
Dilution Range	Peptide ELISA: antibody detection limit dilution 1:128000. IHC: Paraffin embedded Human Spleen. Recommended concentration: 5µg/ml. IF: Strong expression of the protein seen in the nuclei and cytoplasm of HepG2 and U2OS cells. Recommended conc
Formulation	0.5 mg/ml in Tris saline, 0.02% sodium azide, pH7.3 with 0.5% bovine serum albumin. NA
Isotype	IgG
Storage Instruction	Store at -20°C on receipt and minimise freeze-thaw cycles.

TARGET INFORMATION


Gene ID	7454
Gene Symbol	WAS
UniProt ID	WASP_HUMAN
Immunogen Region	Internal
Immunogen Sequence	SPADKKRSGKKKI
Specificity	No cross-reactivity expected with N WASP (WASL).

ADDITIONAL INFORMATION


Note STRICTLY FOR FURTHER SCIENTIFIC RESEARCH USE ONLY (RUO). MUST NOT TO BE USED IN DIAGNOSTIC OR THERAPEUTIC APPLICATIONS.

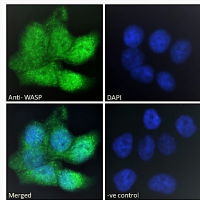

STJ70672 Immunofluorescence analysis of paraformaldehyde fixed HepG2 cells, permeabilized with 0. 15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml) , showing nuclear and cytoplasmic staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

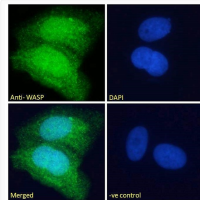

STJ70672 Immunofluorescence analysis of paraformaldehyde fixed U2OS cells, permeabilized with 0. 15% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml) , showing nuclear and cytoplasmic staining. The nuclear stain is DAPI (blue). Negative control: Unimmunized goat IgG (10ug/ml) followed by Alexa Fluor 488 secondary antibody (2ug/ml).

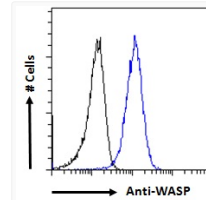

STJ70672 Flow cytometric analysis of paraformaldehyde fixed HepG2 cells (blue line) , permeabilized with 0. 5% Triton. Primary incubation 1hr (10ug/ml) followed by Alexa Fluor 488 secondary antibody (1ug/ml). IgG control: Unimmunized goat IgG (black line) followed by Alexa Fluor 488 secondary antibody.